## NOTICE D'INSTRUCTION - GANT DE BOITE A GANT EN POLYURETHANE (PUR) 6/10mm MONTE SUR BAGUE SUPPORT

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## INSTRUCTIONS - GLOVEBOX GLOVES MADE OF POLYURETHANE (PUR) 6/10mm MOUNTED ON SUPPORT RING

➤ **Scope/Characteristics of the material and expiry**

- This glove can be used in nuclear industries, pharmaceutical industries and miscellaneous industries.
- This glove is manufactured using polyurethane
- This glove protects against certain mechanical risks and against radioactive contamination.
- It can be used in an atmosphere containing ozone. This glove does not protect against ionising radiation.
- No protection against chemical risks is claimed.
- The shelf life of this glove in its original packaging stored under the conditions described below is 4 years.

➤ **Legislation, Risk analysis and recommendations for use**

- This glove complies with regulation 2016/425
- Please note, according to the regulation (UE) 2019/1832, the necessity to carry out a risk analysis relating to wearing the glove (e.g. breaking of tightness)
- Before each use, it is the responsibility of the user to check the quality and the integrity of the glove. In case of any doubt, it is recommended to discard the gloves
- For the assembly and disassembly of the gloves on the gloveboxes, please follow the instructions drafted by the safety manager and the glove port manufacturer.
- In case of a contamination, follow the instructions drafted by the safety manager.
- The gloves should be used at ambient temperature. Please contact the manufacturer for use under other temperature conditions.
- The gloves must be worn on clean and dry hands, with the nails cut short. Avoid wearing jewellery.
- It is recommended to not wear gloves when there is a risk of being caught by moving machine parts.
- This glove can only be mounted if the client has a J2L system that allows receiving the glove (System not validated by IFTH during the EC type examination)

➤ **Components:**

- This glove does not contain substances in proportions in which they are known or suspected to have harmful effects on the user's health or hygiene under the foreseeable conditions of use.

➤ **Properties of the glove**

- **Mechanical Strength EN388+A1 (2018)**

Tests	Level	Pictogram
Abrasion	4 of 4	 4122X
Cut by slicing	1 of 5	
Tear	2 of 4	
Perforation	2 of 4	
Tests Section as per EN ISO 13997 :1999	x	
X: test not carried out		

- Precision, Integrity and ozone cracking

Tests	Level
Precision <b>EN420+ A1:2009</b>	5 of 5
Integrity (pressure 30 mbar) <b>EN421-2010</b>	Compliant
Ozone cracking resistance <b>EN421-2010</b>	4 of 4

The highest class corresponds to the highest performance level

# INSTRUCTIONS - GLOVEBOX GLOVES MADE OF POLYURETHANE (PUR) 6/10mm MOUNTED ON SUPPORT RING

➤ **Marking**

- Example of marking put on the glove and the packaging (see below)

**PIERCAN LOGO**

**PIERCAN (1)**


**0333 (2) 57800 PUR 10 E6 GJABS [XXXXXX]**

(3) (4) (5) (6) (7) (8) (9)

**FR XXXXXX BATCH:XX XX (10)**

**EXPIRY : XX/XX/XX (11)**

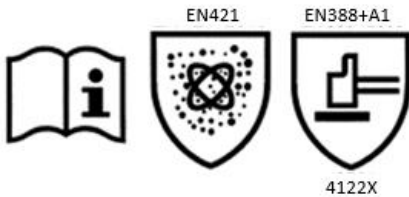(12)

(13)

(14)

No.	DESCRIPTION
1	Manufacturer – In charge of release to the market
2	“CE” marking and No. of the body in charge of the annual inspection
3	Form
4	Length in mm
5	Material
6	Size
7	Thickness in 10th of mm
8	Glove mounted on ring
9	Product Item Code
10	Case number and batch number (“FR” French manufacturing “US” American manufacturing)
11	Date of expiry in storage
12	“INFORMATION” pictogram
13	“Protection against radioactive contamination” pictogram EN 421-2010
14	“Protection against mechanical risks” pictogram EN 388+A1 (2018)

➤ **Size available in the main ambidextrous model**

FORM	Ø RDG (glove port diameter)	Ø BS (internal diameter of the support ring)	SIZE					
			7.5	8.5	9	9.5	10	10.5
64 (GJABS)	136	131					x	

➤ **Suitable packaging for transport**

- The gloves should be transported in their original packaging.

➤ **Storage**

- Store in a dry place, away from light
- Storage temperature: 5 to 25°C.
- Do not store near electrical installations to avoid accelerated ageing

➤ **Cleaning/decontamination and maintenance**

- If necessary, the gloves are washed with soapy water, rinsed with clear water, and dried under a low temperature air flow (<40°C) and powdered. The performance is not guaranteed by the IFTH after washing
- Do not use chemical products
- If it is necessary to use chemical products, please contact the manufacturer for more information.
- Do not use pointed or sharp objects like wire brushes, sandpaper or similar objects

➤ **Treatment of the glove**

- Non-recyclable gloves
- It can be discarded as simple non-toxic waste if it is not soiled by a hazardous product
- If a glove is soiled, discard it through the appropriate channel

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