

➤ Scope/Characteristics of the material and expiry

- This glove can be used in the medical/pharmaceutical industries (domain of life sciences) and nuclear industries.
- This glove is manufactured using natural rubber (polyisoprene).
- This glove protects against certain chemical and mechanical risks, against radioactive contamination and provides protection against micro-organisms and viruses.
- It can be used in an atmosphere containing ozone. This glove does not protect against ionising radiation.
- The shelf life of this glove in its original packaging stored under the conditions described below is 2 years.

➤ Legislation, Risk analysis and recommendations for use


- This glove conforms to regulation 2016/425
- Please note, according to directive (UE) 89/656/EEC, the necessity to carry out a risk analysis relating to wearing the glove (e.g. breaking of leak tightness)
- Claimed performances are guaranteed only if the glove is fitted respecting the direction of delivery of the glove in the bag.
- Before each use, it is the responsibility of the user to check the quality and the integrity of the glove. If the glove shows any tears, holes or changes in surface appearance or color that may indicate tampering with chemicals, the glove should be discarded.
- For the assembly and disassembly of the gloves on the gloveboxes, please follow the instructions drafted by the safety manager and the glove port manufacturer.
- In case of a contamination, follow the instructions drafted by the safety manager.
- The gloves should be used at ambient temperature. Please contact the manufacturer for use under other temperature conditions.
- The gloves must be worn on clean and dry hands, with the nails cut short. Avoid wearing jewellery
- In order to prevent possible risks of allergy and guarantee hand hygiene, it is recommended to wear disposable pre-gloves and to cover the forearm (long sleeves, cuffs, etc.)
- It is recommended to not wear gloves when there is a risk of being caught by moving machine parts.
- During the use of unspecified chemical products (those not mentioned in the list of chemical products), please contact the manufacturer for more information.
- This glove can only be mounted if the client has a J2L system that allows receiving the glove (System not validated by IFTH during the EC type examination).
- The Leakage is not guaranteed by PIERCAN if the ring and glove are separated

➤ Components:

- This glove is likely to induce allergies in sensitive persons owing to the presence of natural rubber proteins. However, it does not contain other substances in proportions in which they are known or suspected to have harmful effects on the user's health or hygiene under the foreseeable conditions of use.


➤ Properties of the glove

- *Protection against chemical products evaluated as per EN ISO 374-1+ A1 :2018*
- *Resistance to deterioration evaluated as per EN374-4 :2013*

Tested products		Level of permeation (palm and cuff)	Average deterioration –palm (%)	Average deterioration – cuff (%)	Pictogram
A	Methanol	1 of 6	18.6	19.5	
N	Acetic acid 99%	2 of 6	21.0	18.3	
P	Hydrogen peroxide 30%	6 of 6	-23.7	-41.6	

This information does not reflect the actual duration of protection at the workplace, or the differentiation between the mixtures and pure chemical products. The chemical resistance has been evaluated under laboratory conditions using the samples collected in the palm and the cuff and only concerns the tested chemical product. It can be different if it is used in a mixture. It is recommended to check that the gloves are suitable for the intended use because the conditions at the place of work can differ from the typical test conditions, depending on the temperature, abrasion and deterioration. When they are worn, the protective gloves provide less resistance to hazardous chemical products owing to the change in their physical properties. Movements, tears, friction or deterioration caused by contact with chemical products, etc. can considerably reduce the actual duration of use. For corrosive chemical products, deterioration can be the most important factor to be taken into account in the selection of chemical-resistant gloves. Before use, it is recommended to inspect the gloves in order to ensure that they have no defects or imperfections.

- *Protection against viruses, bacteria and moulds evaluated as per EN ISO 374-5: 2016*
- *Leak-tightness evaluated as per EN374-2 :2014*

Tests	Level	Pictogram
Airtightness	Compliant	 VIRUS
Water-tightness	Compliant	
Penetration by contaminated liquids under hydrostatic pressure (ISO 16604 method B)	No passage	

The resistance to penetration has been evaluated under laboratory conditions and only pertain to the tested test specimen.

Attestations d'Examen UE de Type (AET) (EU type examination certifications) delivered by: I.F.T.H. (no. 0072) Avenue Guy de Collongue – 69134 ECULLY Cedex.

Monitoring organisation: AFNOR CERTIFICATION (no. 0333):11, rue Francis Pressensé-93571 LA PLAINE ST DENIS Cedex


PIERCAN: Impasse des Macareux ZI Huppain / 14520 PORT EN BESSIN FRANCE Tel.: 33 (0)2 31 21 73 80 Fax: 33 (0) 2 31 21 40 23 Email: piercan @ piercan.fr Website: www. Piercan.fr

Instructions: Edition of 30/06/2025 Rev 4



INSTRUCTIONS - GLOVEBOX GLOVES MADE OF NATURAL RUBBER (H) 6/10mm MOUNTED ON SUPPORT RING

- Protection against mechanical risks evaluated as per **EN388+A1: 2018**

Tests	Level	Pictogram
Abrasion	X of 4	 XX11X
Cut by slicing	X of 5	
Tear	1 of 4	
Perforation	1 of 4	
Tests Section as per EN ISO 13997 :1999	X	
X: test not carried out		

- Precision, Integrity and ozone cracking

Tests	Level
Precision EN ISO 21420 + A1 :2024	5 of 5
Integrity (pressure 30 mbar) EN421-2010	Compliant
Ozone cracking resistance EN421-2010	1 of 4

The highest figures correspond to the highest performance levels.

➤ Integrity EN421-2010

The “glove port diameter” and “form” lists are not exhaustive. For any requests, please contact PIERCAN.

GLOVE PORT DIAMETER PORT DIAMETER (mm)	FORM – MODEL (mm)	ORIGINAL LEAK-TIGHTNESS PRESSURE AIR LEAK TEST INITIAL PRESSURE (mbar)
136	14	30
156	10	30

➤ Marking

- Example of marking put on the glove and the packaging (see below)

LOGO PIERCAN




PIERCAN

(1)


10750 H 9 E6 GJABS [XXXXX]

(3) (4) (5)(6)(7) (8) (9)

FR XXXXXX LOT : XXXX  AAAA/MM

(10)

(11)

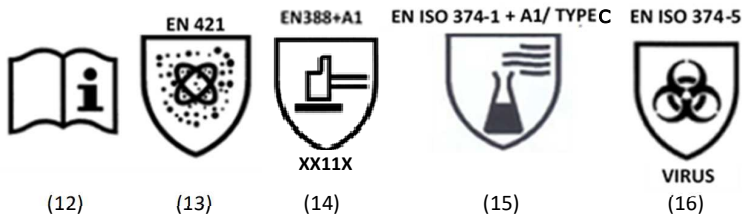

(12)

(13)

(14)

(15)

(16)

No.	DESCRIPTION
1	Manufacturer – In charge of release to the market
2	“CE” marking and No. of the body in charge of the annual inspection
3	Form
4	Length in mm
5	Material
6	Size
7	Thickness in 10th of mm
8	Glove mounted on support ring
9	Product Item Code
10	Case number and batch number (“FR” French manufacturing “US” American manufacturing)
11	Date of shelf life in storage + pictogram
12	“INFORMATION” pictogram
13	“Protection against radioactive contamination” pictogram EN 421 :2010
14	“Protection against mechanical risks” pictogram EN 388+A1:2018
15	“Protection against chemical risks” pictogram EN ISO 374-1 +A1 :2018
16	“Protection against micro-organisms and viruses” pictogram EN ISO 374-5 :2016

➤ Size available in the main ambidextrous model

FORMS	Ø RDG (Glove port diameter)	Ø BS (Internal diameter of the support ring)	SIZES				
			6	7	8	9	10
14 (GJA+BS)	136	131		x	x	x	x
10 (GJA + BS)	156	156		x	x	x	

Attestations d'Examen UE de Type (AET) (EU type examination certifications) delivered by: I.F.T.H. (no. 0072) Avenue Guy de Collongue – 69134 ECULLY Cedex.

Monitoring organisation: AFNOR CERTIFICATION (no. 0333):11, rue Francis Pressensé-93571 LA PLAINE ST DENIS Cedex

PIERCAN: Impasse des Macareux ZI Huppain / 14520 PORT EN BESSIN FRANCE Tel.: 33 (0)2 31 21 73 80 Fax: 33 (0) 2 31 21 40 23 Email: piercan @ piercan.fr Website: www. Piercan.fr

Instructions: Edition of 30/06/2025 Rev 4



INSTRUCTIONS - GLOVEBOX GLOVES

MADE OF NATURAL RUBBER (H) 6/10mm MOUNTED ON SUPPORT RING

➤ **Suitable packaging for transport**

- The gloves should be transported in their original packaging.

➤ **Storage**

It is recommended to store gloves:

- In their original packaging, flat, black side of bag up and in their original box
- At a storage temperature between 5 to 35°C /41 to 95°F
- In a dry place away from direct light
- Away from electrical power source to prevent accelerated ageing.

➤ **Cleaning/decontamination and maintenance**

- Gloves are not designed to be laundered.
- Prior to the use of unspecified chemicals, please contact the manufacturer for compatibility information.
- Do not contact gloves with sharp or pointed objects such as wire brushes, sandpaper, screwdrivers, or similar objects.
- If contaminated with chemicals, the gloves are for single use only

➤ **Treatment of the glove**

- Non-recyclable glove
- It can be discarded as simple non-toxic waste if it is not soiled by a hazardous product
- If a glove is soiled, discard it through the appropriate channel

The EU Declaration of Conformity is available at: www.piercan.fr